


Introduction to X-Ray Image Enhancement System

❏ Made by Nivethitha Chowthri, Monisha D , Nithisha Paulin

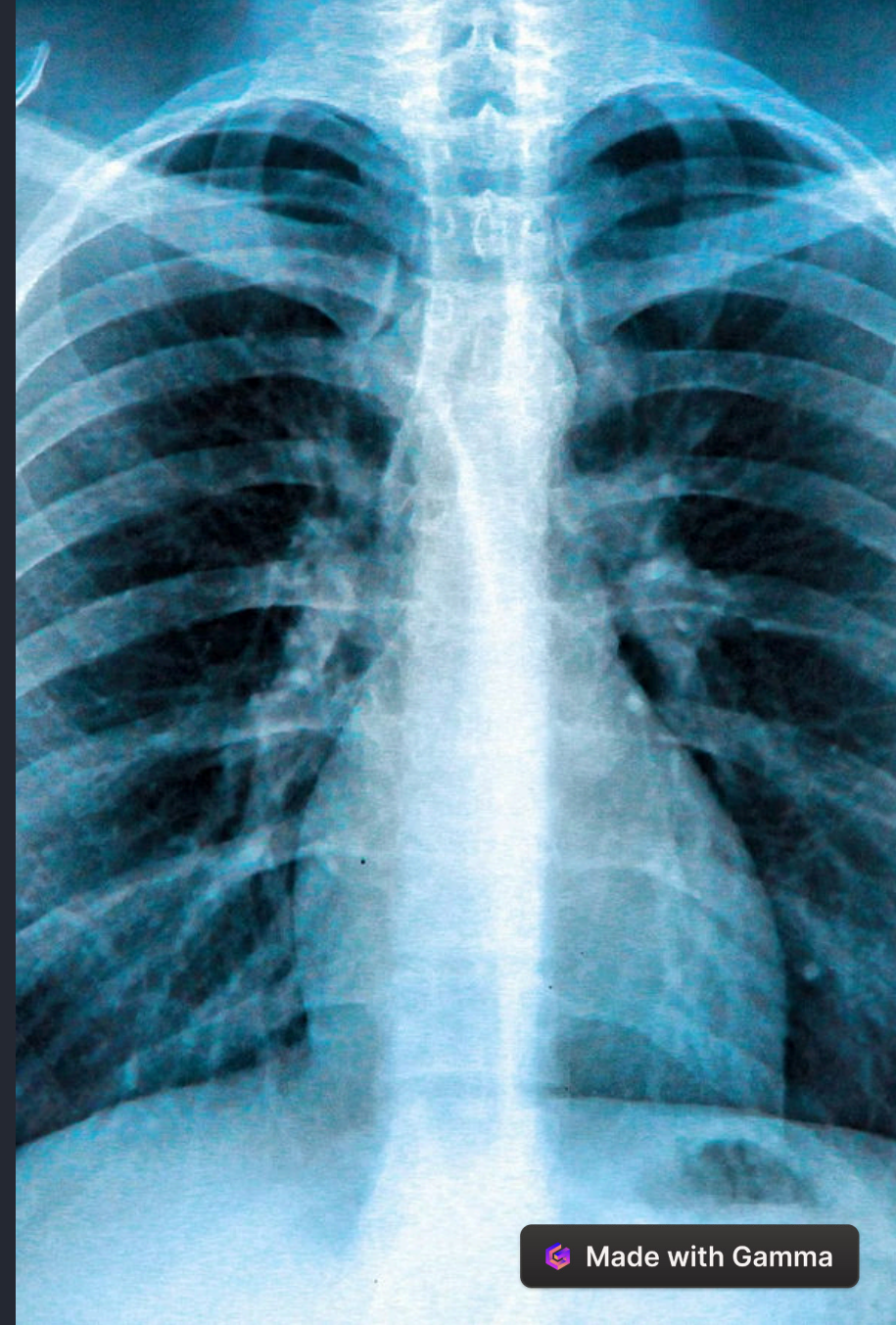

Challenges in X-Ray Image Quality

1 Low Contrast


X-ray images can suffer from low contrast, making it difficult to distinguish between different tissues and structures.

2 Noise and Artifacts

X-ray images are susceptible to various types of noise and artifacts that can obscure important details.

3 Varying Tissue Densities

The wide range of tissue densities captured in a single x-ray image can make it challenging to optimize visibility for all regions.

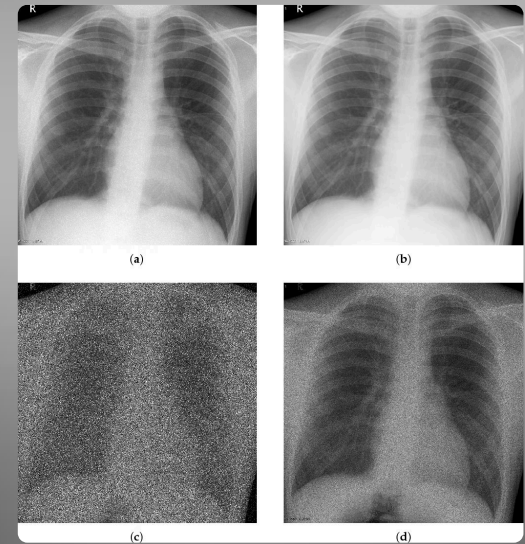

Some of the X-ray images :



UM VS HEF VS CLAHE

Unsharp Masking (UM) :

Unsharp masking is one of a technique used in image processing to improve the visual clarity and highlight important details. Unsharp masking is a linear filter that is capable of amplify high frequencies of an image. For blurring the image it uses Gaussian Filter, Median Filter, Maximum Filter, Minimum Filter.

Steps Involved in UM:

1. Creates a Blurred Version of the Original Image
2. Subtract the Blurred Image from the Original Image
3. Add the Difference Back to the Original Image

$$\text{sharpenedimage} = \text{originalimage} + \text{amount} * \text{unsharpedmask}$$

$$\text{unsharpedmask} = \text{originalimage} - \text{blurredimage}$$

Results:



High-frequency Emphasis filtering (HEF) :


High-frequency Emphasis filtering is a technique that uses Gaussian High Pass Filter to emphasize and accentuate the edges. The edges tend to be expressed in the high-frequency spectrum since they have more drastic changes of intensity.

Steps involved:

1. It applies a gaussian high pass filter into it.
2. The image has to go through the Fourier transformation and the filter function is calculated onto it.
3. After the inverse transformation we will have filtered image. Secondly, the contrast of the image will be adjusted with simple Histogram Equalization

$$\text{sharpened image} = (\text{original image} + (\text{Gaussian Highpass Filter})) * (\text{Histogram Equalization})$$

Results of HEF:

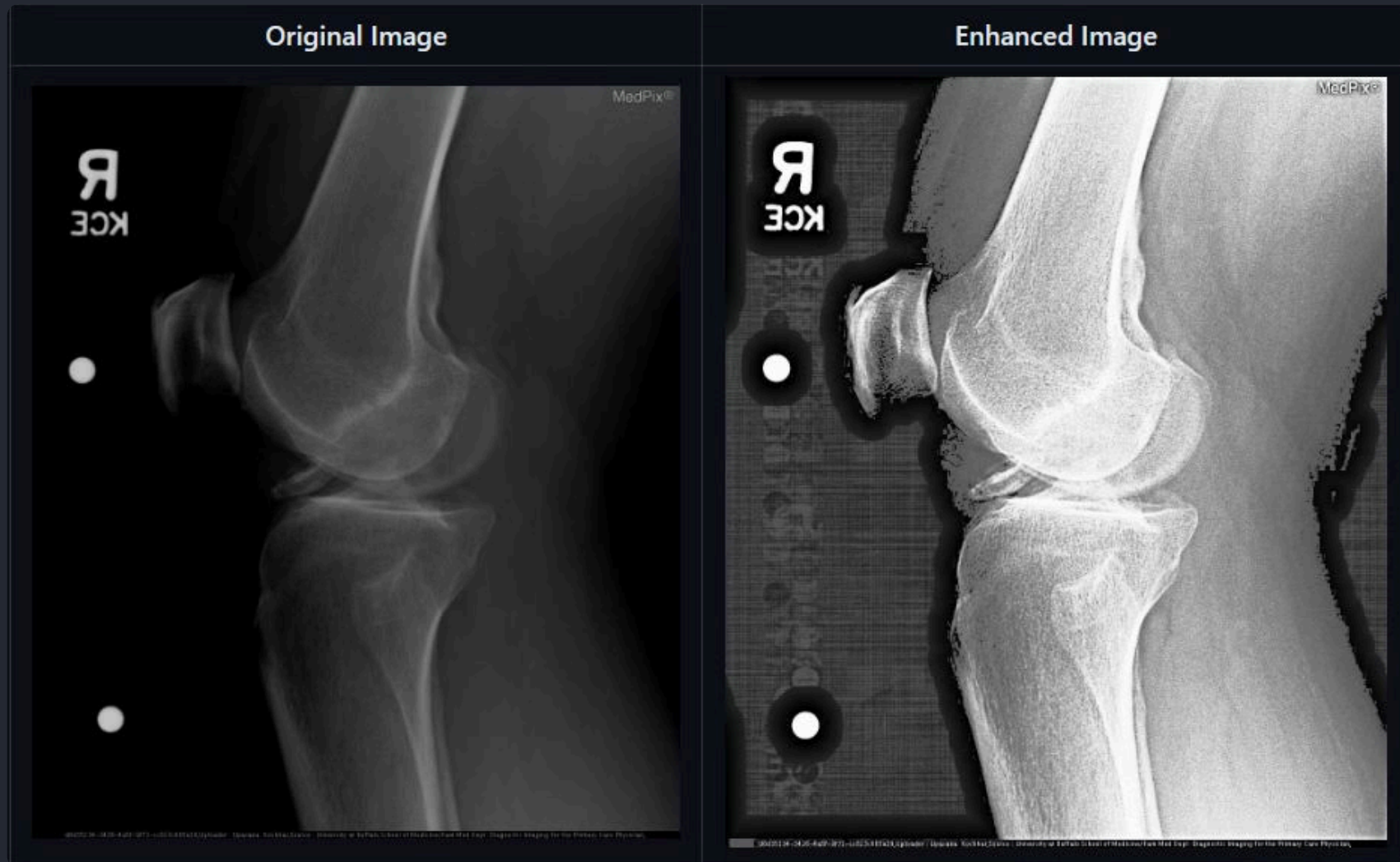

Contrast Limited Adaptive Histogram Equalization (CLAHE):

Contrast Limited Adaptive Histogram Equalization is the most advanced image enhancement technique used to improve the contrast of images. It enhances the local contrast of an image, making details more visible, while preventing over-amplification of noise.

Steps Involved:

1. Convert the images to grayscale and then normalize
2. Divide the Image into Small Regions
3. Apply Histogram Equalization to Each Tile
4. Clip the Histogram
5. Redistribution of Clipped Pixels
6. Bilinear Interpolation
7. Combine the Tiles

Implementation of CLAHE for X-Ray Images

Window Size Optimization

The size of the rectangular region around the pixel to be processed

Clip Limit Adjustment

The clip limit parameter controls the amount of contrast enhancement, and it should be tuned based on the specific characteristics of the x-ray image.


Iterative Processing

Applying CLAHE iteratively can further improve the image quality by gradually enhancing the contrast and reducing noise.

Image-Specific Tuning

The CLAHE algorithm may require different parameter settings for different types of x-ray images to achieve optimal results.

Results of CLAHE:

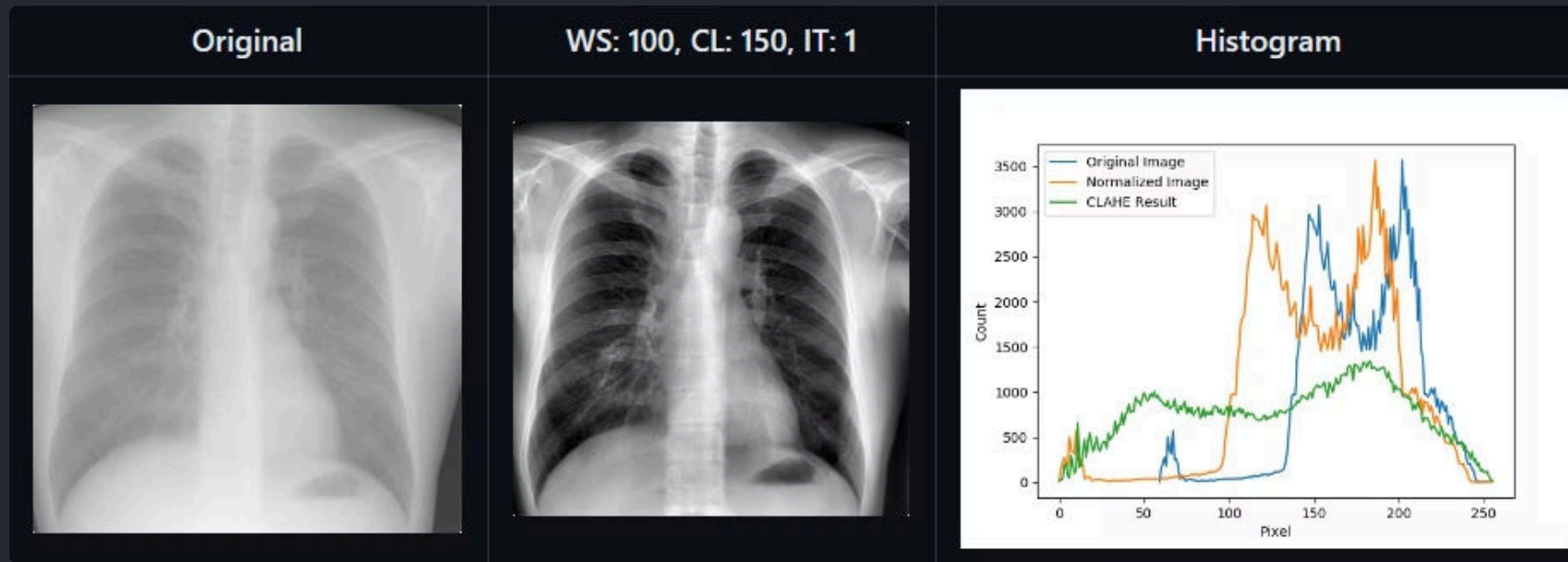

Advantages of Using CLAHE for X-Ray Enhancement



Improved Contrast

CLAHE effectively enhances the contrast of x-ray images, making it easier to visualize and interpret the underlying structures.



Noise Reduction

The contrast limiting feature of CLAHE helps to suppress noise and artifacts, resulting in cleaner and more informative x-ray images.



Adaptability

CLAHE can be tailored to work well for a variety of x-ray imaging modalities and anatomical regions by adjusting the algorithm parameters.




Computational Efficiency

The CLAHE algorithm is relatively fast and can be implemented in real-time, making it suitable for clinical applications.

Comparison to Other Enhancement Techniques

Technique	Contrast Improvement	Noise Suppression	Computational Complexity	Run time
HEF	Low	High	Medium	Medium Fast
Unsharp Masking	Medium	Medium	Medium	Fastest
CLAHE	High	High	High	Slowest

Comparison of UM, HEF, CLAHE (images)

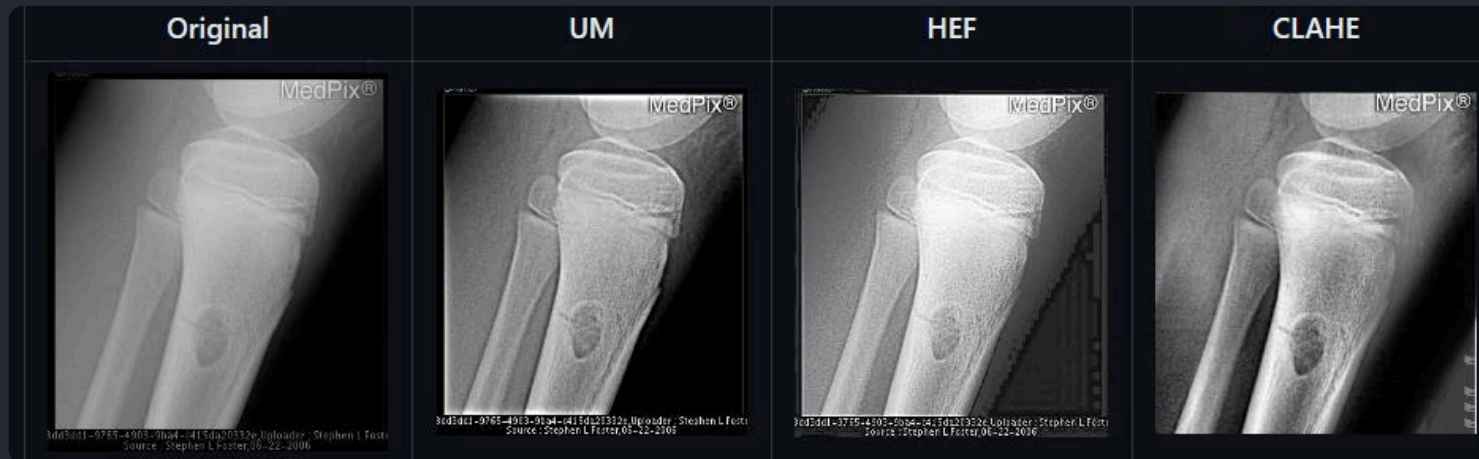

Conclusion and Future Enhancement

1

Improved Diagnostics

CLAHE-enhanced x-ray images provide better visibility and clarity, leading to more accurate diagnoses and improved patient outcomes.

2

Automated Analysis

The enhanced image quality from CLAHE can enable more sophisticated computer-aided detection and analysis algorithms.

3

Personalized Optimization

Future research may focus on developing adaptive CLAHE algorithms that can automatically optimize parameters for each individual patient's needs.

4

Improved Runtime

With the help of artificial intelligence, we can improve the run time even faster than the previous model.



(b)



(e)



(h)

THANK YOU